Fiber-Lite[®]

Dolan-Jenner[®]

Model Mi-LED LED Fiber Optic Illuminator

Dolan-Jenner introduces the New Fiber-Lite® Mi-LED, LED Fiber Optic Illuminator for microscopy and general illumination. The Mi-LED utilizes the latest technology and is energy efficient, featuring cool white LED lighting with a "customer requested" color temperature of 5000K and even higher output than the Fiber-Lite Mi-150, 150 watt halogen fiber optic illuminator.

Microscopy Applications


The Fiber-Lite[®] Mi-LED was designed with microscopy applications in mind, able to support and fully illuminate gooseneck and annular ring light fiber optics. It is also a great solution for inspection equipment, forensics, machine vision, industrial borescopes, assembly stations, and OEM custom lighting.

LED Lighting

The LED lighting of the Fiber-Lite® Mi-LED is an energy efficient cost saving alternative requiring less power while delivering higher performance than comparable 150 watt halogens. Offering even more value, the Fiber-Lite® Mi-LED is practically maintenance free with no lamps to fail, and no lamp sockets to wear out. The long life LED's of Fiber-Lite® Mi-LED are rated to provide over 40,000 hours of light. Conversely, a 150 watt halogen can only provide around 200 hours of light. With such a dramatic increase in light hours, large savings can be realized through greater up time during applications, lower maintenance, and the elimination of replacement bulbs.

Accessories

The B1 (15mm), A2 (25mm) and A7 (iris) nosepiece adapters make this a drop in replacement and compatible with all standard Dolan-Jenner Mi-150 and 180 series fiber optic light guide configurations. Furthermore, the Fiber-Lite® Mi-LED has a universal voltage input of 100-240 VAC, and is RoHS compliant & CE Certified.

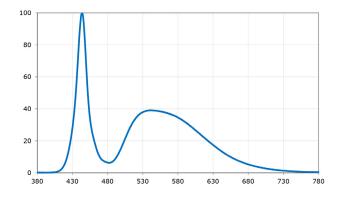



"WARNING: This product can expose you to carbon black, acrylonitrile, or ethylbenzene which are known to the State of California to cause cancer (For more information go to www.p65warnings.ca.gov)"

- High Power White LEDs
- Consistent Color Temperature
- Made in the U.S.A.

Model Mi-LED Features:

- Universal Power Supply
- Low (IR) Heat
- Low Energy Consumption
- Stackable Housing
- Multiple Fiber Optic Interface Options
- Small Footprint
- 2 Year Warranty

Configurations:

- 15mm ID accepts "Mi-150" light guides
- 25mm ID accepts "SX" and "MX" series
- 25mm ID with manual Iris accepts "SX" series adapters

Model Mi-LED LED Fiber Optic Illuminator

Dolan-Jenner[®]

ORDERING INFORMATION

MI-LED-	_	-		
Model	Corc	l	0pti	on
MI-LED	US	North American Cord	B1	15mm
	EU	European Cord 230V	A2	25mm
	UK	United Kingdom 230V	A7	25mm w/ Iris
			DG	Dual Gooseneck Fiber
			SG	Single Gooseneck Fiber
			RA	Nikon SMZ Ring Light
			RB	66mm Ring Light
			RD	Olympus SZ Ring Light
			RG	Leica S4 Ring Light
			RL	60mm Ring Light

GENERAL SPECIFICATIONS

Physical Description		Performance Data		
Dimension	7.0"W x 4.9"H x 8.8"DP	Light Output	2000+ Lumens (12.7mm aperture)	
Weight	6 lbs.	Lighting	High Power White LEDs	
Intensity Control	Variable 10 - 100%	Noise Rating	<30 dB	
Lighting Properties		Electrical Data		
LED Life 40,000 Hours D70		Input Voltage	100-240 VAC 50/60 Hz	
Color Temp. Typical: 5000 K		Environmental Data		
Certifications		Operating Temp.	5 - 40°C	
ETL, CE, RoHS		Cooling	Fan Cooled	

OPTIONAL FIBER OPTIC INTERFACES



SYSTEMS

